CAT CYSTATIN-C SPARCL™ ASSAY Life Diagnostics, Inc., Catalog Number: CYSC-SP-8

INTRODUCTION

Cystatin-C, a protein of ~13 kDa, is used as a biomarker of kidney function in humans. Recent studies indicate that it may serve a similar purpose in cats (refs 1-2).

PRINCIPLE OF THE ASSAY

The cat cystatin-C SPARCL™1 (Spatial Proximity Analyte Reagent Capture Luminescence, ref 3) assay uses two cat cystatin-C specific antibodies. One is conjugated to horseradish peroxidase (HRP), the other is conjugated to acridan, a chemiluminescent substrate. When the HRP and acridan conjugated antibodies bind to cystatin-C they are brought into proximity. With the addition of hydrogen peroxide, HRP catalyzes oxidation of proximal acridan molecules causing a flash of chemiluminescence. Acridan conjugated antibodies distant from HRP produce no signal. This principle allows the development of a homogeneous assay that allows rapid measurement of cystatin-C concentrations.

Diluted serum samples and standards are mixed with the HRP and acridan-conjugated antibodies in the wells of the 96-well SPARCL™ plate² provided with the kit. After incubation for 30 minutes on a shaker at 25°C and 150 rpm, the plate is placed into a luminometer. Trigger solution containing hydrogen peroxide is injected into each well and luminescence is immediately measured. The concentration of cystatin-C is proportional to luminescence and is derived from a standard curve.

MATERIALS AND COMPONENTS

Materials provided with the kit:

Anti-cystatin-C HRP stock

Store ≤ -70°C

Anti-cystatin-C acridan stock

Store ≤ -70°C

Cystatin-C stock.

Store ≤ -70°C

- Diluent; CSD50-1, 2 x 50 ml
- Trigger solution; TS7-1, 7 ml
- White SPARCL™ plate (12 x 8-well)
- Clear untreated 96-well plate

Materials required but not provided:

- Precision pipettes and tips
- Polypropylene microcentrifuge tubes
- Vortex mixer
- Plate incubator/shaker
- Luminometer capable of simultaneous injection/measurement
- Curve fitting software

STORAGE

Store the HRP conjugate, acridan conjugate and cystatin-C stock at -70°C (they may be stored at -20°C for one week). The remainder of the kit should be stored at 2-8°C. The SPARCL™ plate should be kept in a sealed bag with desiccant and antioxidant. The kit will remain stable for at least six months from the date of purchase, provided that the components are stored as described.

GENERAL INSTRUCTIONS

- Please take the time to completely read all instructions before starting your assay. Contact us if you need clarification.
- All reagents used in the assay should be allowed to reach room temperature (25°C) before use.
- It is important that standards and samples be added to the SPARCL[™] plate quickly. If testing many samples, rather than pipetting standards and samples directly into the white SPARCL™ plate using a single channel pipettor, we recommend the following. First, pipette an excess volume of standards and samples into appropriate wells of the clear 96well plate. Then use an 8- or 12-channel multipipettor to guickly and efficiently transfer 50 µl aliquots to the appropriate wells of the white SPARCL™ plate. The wells of the clear plate have a maximum volume of 300 µl.
- Follow the sequence of events below when running the assay.

Prime and program the Luminometer

Prepare standards and diluted samples

Û

Prepare HRP + Acridan conjugate mix

①

Add HRP + Acridan conjugate mix to the wells (25 μ l) Û

Add standards and samples to the wells (50 μ l)

Û

Incubate plate at 150 rpm/25°C for 30 min

Ú

Measure luminescence after injection of Trigger (37.5 μl)

STANDARD PREPARATION

The cystatin-C stock is comprised of cat serum of known cystatin-C concentration diluted in a stabilizing carrier protein matrix³. Thaw the stock shortly before use.

- Label five polypropylene tubes as 1, 0.5, 0.25, 0.125 and 0.0625
- Into the tube labeled 1 ng/ml, pipette the volume of diluent 2. detailed on the cystatin-C stock vial label. Then add the indicated volume of cystatin-C stock and mix gently. This provides the 1 ng/ml standard.
- Dispense 150 µl of diluent into the tubes labeled 0.5, 0.25, 0.125 and 0.0625 ng/ml.
- Pipette 150 µl of the 1 ng/ml cystatin-C standard into the tube labeled 0.5 ng/ml and mix. This provides the 0.5 ng/ml cystatin-C standard.
- Similarly prepare the remaining standards by serial dilution.

Please Note: Use the standards within one hour of preparation.

SAMPLE PREPARATION

Serum samples should be diluted at least 50-fold with diluent to avoid matrix effects. Plasma and hemolyzed serum should not be used. In

¹ The SPARCL technology was developed by Lumigen Corp.

² The plate provided with the kit has been treated with a reagent that reduces background chemiluminescence. Untreated plates cannot be used.

³ The concentration of cystatin-C was determined using purified human cystatin-C as reference standard.

studies at Life Diagnostics, we found serum levels of 230.0±89 ng/ml (mean±SD, n=6). Optimal dilutions must be determined by the end user; however, we suggest an initial dilution of 400-fold (1.25 μl of serum mixed with 498.75 μl of diluent). Use only the diluent provided with the kit.

CONJUGATE MIX PREPARATION

Instructions for preparation of the conjugate mix are detailed on the box that contains the HRP and acridan conjugates. If necessary, after thawing, briefly centrifuge to ensure that the contents are at the bottom of the tubes. Prepare the mix shortly before use using the diluent provided with the kit.

LUMINOMETER SETUP

- 1. The luminometer must be capable of injection and simultaneous measurement of luminescence without any delay.
- 2. Prime the luminometer injection port with 1 ml of trigger solution.
- Place the injection needle into the injection port as needed for BMG luminometers.
- 4. Program the luminometer to inject 37.5 μ l of trigger solution per well and to measure from time zero for 1 second (50 x 0.02 second intervals).
- 5. Define the format of the assay using the luminometer software.
- 6. Because the white SPARCL™ plate is provided as a 12 x 8-well strips, allowing use of fewer than 96-wells, make sure that the luminometer is programmed to inject trigger solution only into the wells being used.
- 7. We use a BMG LUMIstar Omega set at a gain of 3600. Optimal gain should be determined by the end user.
- There are a number of manufacturers of luminometers that are equipped to run a SPARCL™ assay. Please contact Life Diagnostics or Lumigen (www.lumigen.com) to discuss your luminometer.

PROCEDURE

- Before starting the assay ensure that the luminometer is primed with trigger solution and that the injection needle is positioned in the injection port.
- 2. Secure the desired number of SPARCL™ 8-well strips in the holder. Immediately seal unused strips in the resealable bag with desiccant and antioxidant. Store unused strips at 2-8°C.
- 3. Aliquot 25.0 µl of conjugate mix into each well.
- 4. Dispense $50.0 \,\mu l$ of standards and diluted samples into the wells (we recommend that standards and samples be tested in duplicate).
- 5. Incubate on an orbital micro-plate shaker at 150 rpm and 25°C for 30 minutes.
- After the 30-minute incubation, place the plate in the luminometer and measure luminescence after injection of trigger solution (37.5 μl).
- 7. Remove the plate from the luminometer and discard the used strips. Keep the plate frame if future use is intended.

CALCULATION OF RESULTS


- Before calculating results, review the raw data. If artefacts (RLU spikes) are apparent immediately after injection of trigger solution, eliminate that portion of the luminescence profile from analysis for all wells. We routinely use the sum of RLU values from a 100-980 ms data collection window.
- 2. Determine the sum of RLU values within the data collection window for the standards and samples.

- 3. Using graphing software, construct a standard curve by plotting the sum of the RLU values for the standards versus the log₁₀ of the cystatin-C concentration and fit to a sigmoidal, 4PL model.
- 4. Derive the corresponding concentration of cystatin-C in the samples from the standard curve (remember to derive the concentration from the antilog).
- 5. Multiply the derived concentration by the dilution factor to determine the concentration of cystatin-C in the original sample.
- If the sum of the RLU values of diluted samples fall outside the standard curve, samples should be appropriately diluted and retested.

TYPICAL STANDARD CURVE

A typical standard curve is shown below. This curve is for illustration only and should not be used to calculate unknowns.

Cystatin-C	RLU
(ng/ml)	
1	109782
0.5	52075
0.25	22323
0.125	7017
0.0625	3608

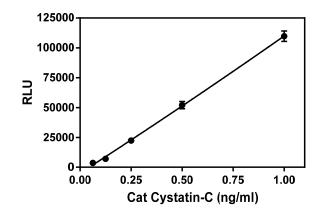

REFERENCES

- Poswiatowska-Kaszczyszyn I. Usefulness of serum cystatin-C measurement for assessing renal function in cats. Bull Vet Inst Pulawy. 56:235-239 (2012)
- 2. Kovarikova S. Urinary biomarkers of renal function in dogs and cats: A review. Veterinarni Medicina. 60:589-602 (2015)
- Akhavan-Tafti H. et al. A homogeneous chemiluminescent immunoassay method. J Am Chem Soc. 20;135(11):4191-4 (2013)

Rev 022218
For technical assistance please email us at techsupport@lifediagnostics.com